



English

Customer Release
Notes

300017891311_A

Advanced Visualization Workspace

Version 16.0.3 2026-08-11

Table of Contents

1	About this document	5
	Disclaimer.....	6
2	Release description	7
3	Main new features and enhancements	11
	V16.0.3	11
4	Resolved issues.....	15
	V16.0.3	15
	AV Viewer	15
5	Known issues.....	17
	V16.0.3	17
	System	17
	Multimodality applications	18
	AV Cardiovascular Suite and AV Viewer	18
	Service & IT	22

1 About this document

This Release Notes document describes new applications and features of Advanced Visualization Workspace V16, AV Viewer V1.0.2, and AV Cardiovascular Suite V1.0.2. This document also lists the resolved issues, known issues, and limitations.

NOTICE

The Release Notes document is not a substitute for the Instructions for Use (IFU) documents. Before using the medical device software, read the IFUs thoroughly, paying particular attention to all warnings, cautions, and notes they contain.

User documentation

Philips' AVW V16.0.3 is accompanied by the following user documentation, including Instructions for Use (IFU) and reference guide documents:

NOTE

Ensure the Instructions for Use correspond with the software version you are operating.

- AVW V16.0.3 System IFU
- AVW V16.0.3 System Administrator Guide
- AV Viewer IFU V1.0.2
- AV Cardiovascular Suite IFU V1.0.2
- AVW V16.0.2 MR Apps IFU
- AVW V16.0.2 Multimodality Apps IFU
- AVW V16.0.3 CT Analysis IFU
- CT ASPECT Scoring V16.0.2 IFU
- NM Suite V16.0.2 IFU
- AVW V16.0.3 Zero Footprint Viewer IFU
- AVW V16.0.3 Spectral Apps IFU
- AVW V16.0.2 CT Review IFU

Instructions for Use and reference guides for third party software launched from AVW16.0.3:

- ADAS 3D IFU
- 3mensio V10.7 IFU
- Caas MR 5.3.3 IFU
- ClearRead V5.8 IFU
- NeuroQuant V5.3 IFU
- NeuroQ V3.8 IFU
- XD 4.2 IFU
- XD 3.6.11-5 IFU

- Cedars-Sinai Cardiac Suite V2017 Rev. I IFU
- ECTb V4.2 IFU
- Corridor4DM V2025 IFU
- XD 3.6.11 Gestures and Shortcuts
- ECTb 4.2 Users Guide
- CS QBS Ref. V2017
- CS AutoQUANT Ref. V2017

NOTE

If there is a version number in the document name, it reflects the latest version that specific document was updated.

NOTE

The list of user documentation above covers the full scope of the offered software options, which are subject to licensing and regulatory approval. As such, some of these documents may not be applicable to you according to your software licenses and the regulatory approvals in your territory.

Disclaimer

Wherever the term **Archive** or its derivatives are used in the **Instructions for Use** or any other document associated with Advanced Visualization Workspace, please treat the meaning as synonymous with the term **Storage**.

2 Release description

The Advanced Visualization Workspace (AVW) is an advanced visualization platform that offers a single integrated solution, spanning clinical domains and modalities, with the power to visualize, diagnose, and communicate with one consistent, guided and automated workflow.

This continuously evolving solution is tailored to the needs of the physicians and the hospital network, designed for seamless integration with other information systems and offers a wide variety of clinical applications.

Advanced Visualization Workspace offers customers the following:

- AVW provides tools to assist in the diagnosis of multiple clinical questions based on multi-modality imaging with capabilities designed for optimized results sharing between the imager and the treating physician.
- AVW offers clinical depth and innovation with comprehensive set of more than 70 clinical applications spanning across clinical domains including oncology, cardiovascular, neurology, pulmonology and others, addressing different clinical needs while focusing on optimizing user experience and efficiency.
- The AVW environment enables launching of CT ASPECTS, NM Suite, AV Viewer, and AV Cardiovascular Suite positioned as SaMD applications.
- AVW offers a solution which is tailored for the enterprise needs and can scale up to fit large scale hospitals and integrated delivery networks (IDNs).
- AVW introduces a single multi-modality and multi-vendor solution which provides a comprehensive single-source view of the patient-related information. AVW allows the patient information to be accessible from anywhere within the organization and also enables maintaining consistent applications, workflow and user preferences.
- Service offering which gives a steady stream of clinical and IT innovation through ongoing software upgrades and updates, allowing to cope with the rapid changes of advanced analysis environment and stay at the forefront of clinical innovation. Clinical training, support on demand and consulting services are offered as part of the end-to-end service offering.
- AVW helps to maximize resources to support productivity through optimized performance, efficient workflows and machine learning capabilities.

Advanced Visualization Workspace V16

Advanced Visualization Workspace version 16 (AVW V16) replaces AVW V15 and introduces advancements in clinical functionality, user experience, interoperability, and workflow automation. AVW V16 introduces AV Viewer and AV Cardiovascular Suite, and integrates enhanced AI capabilities, improved image analysis algorithms, and enterprise interoperability with Philips Image Guided Therapy (IGT) systems. AVW V16 continues Philips' strategy to deliver an integrated platform designed to support precision diagnosis, efficient planning, and collaborative care delivery across modalities and clinical domains.

AVW V16 introduces enhancements to core applications such as:

- CT Liver Analysis
 - AI segmentation of liver and spleen with direct results integration to PACS

- A new AI algorithm for contrast phase classification of multiphase CT
- Segment the hepatic and portal veins on the portal phase.
- Segment the arteries on the early arterial phase.
- Updated workflow steps
- Lung Nodule Analysis
 - Lung-RADS 2022 compliance and manual cyst segmentation and classification tools.
 - For ClearRead CAD users, automatic results are sent to PACS only for images that contain detected findings.
- MR Application Enhancements
 - MR Image Algebra for mathematical operations between source and reference volumes.
 - MR Image Cumulation for generation of new images from multi-echo or b-value sequences.

Cardiovascular Suite

AVW V16 introduces a completely new, all-in-one Cardiovascular Suite designed to address the full spectrum of cardiovascular diagnostic and pre-interventional needs from a single user interface.

Key enhancements include:

- An extensive cardiovascular workflow supports cardiac CT, cardiac MR, and vascular analysis from one access point.
- Cardiac CT Analysis:
 - Automatic coronary segmentation
 - CAD-RADS 2.0-compliant reporting
 - Pre-procedural PCI planning module: a dedicated workflow to support procedural planning
 - Supports CT Spectral data
 - Export results to a defined location
- Cardiac MR:
 - Flexible workflow that supports all cardiac MR analysis tasks
 - Results-driven reporting and customizable layouts per clinical indication
 - Export results to a defined location
- Vascular Analysis:
 - Unified CTA/MRA workflow
 - Aortic AI wall segmentation algorithm for CTA cases
 - New planning tools and template editor for endovascular procedures
 - Export results to a defined location

This suite aims to deliver one intuitive environment for cardiovascular clinicians and supports enterprise scalability and communication across all departments.

AV Viewer

AVW V16 introduces AV Viewer, a single framework for multi-modality visualization and results review.

- Streamlined reading environment across CT, spectral CT, MR, NM, XA and Ultrasound datasets
- Configurable hanging protocols and personalized reading workflows
- Supports basic measurements and side-by-side study comparison
- An integrated Findings Dashboard for consolidated results display and export

New and updated third-party applications

Third-party applications new to AVW V16:

- 3mensio:
 - For structural heart disease
 - Dedicated workflows to assess structures in the heart to plan different procedures

AVW V16 introduces the following third-party application updates:

- New version of CleaRead CAD
- Corridor 4DM V2025:
 - Streamlined amyloidosis quantification with new reporting templates, expanded tracer support, validated iMFR and TPD metrics, and SUV SPECT quantification
 - Optional 4DM Reconstruction module for integrated SPECT image reconstruction and motion correction

Enhancements based on customer feedback

A host of applications have been updated across the portfolio based on customer feedback.

3 Main new features and enhancements

Advanced Visualization Workspace V16 provides new applications, features, and enhancements which are described in the following sections.

V16.0.3

Viewer

- Added the functionality "Add to Running Application", where you can open new studies in an existing viewer session.

Spectral Magic Glass

Interactively explore spectral information within a localized region of interest while maintaining the surrounding anatomical context. The Spectral Magic Glass enables fast, intuitive comparison of spectral results without switching between series, supporting more efficient image review and interpretation. Additional improvements include:

- Save batches of spectral results as original images with burned-in labels.
- New factory hanging protocols for improved workflow consistency.
- Selected spectral curve is now maintained when switching between spectral results.
- SBI pictorial image now displays anatomical reference information.
- SBI series are identified with a dedicated spectral icon for easier recognition.

Control curve visibility in the viewer

Show/Hide the CL curve using a checkbox on the Curve list. When you select a new curve, it is automatically shown.

Grid overlay

Display a grid on Slab and 3D viewports, either using the right-click shortcut menu's Image Options, or the **G** keyboard shortcut. You can adjust grid spacing using the floating Grid Spacing control.

Optimizations to Hanging Protocols

- When loading a study to the viewer, only the valid hangings for the loaded study appear in the dropdown list.
- You can define your favorite Hanging Protocols (HPs) and even define several favorites. For example, for Cardiac MR your favorite hanging may be "Functional", for CT your favorite hanging may be "Slab CT" and for Spectral, a spectral Hanging Protocol.
 - Upon loading a study, the application selects the first valid, favorite HP.
- You can define an HP with fusion, i.e. CT-PET.
- 3D orientation is retained as part of the HP.

Cardiovascular Suite

You can save the motion-corrected source images and map as a batch.

- Cardiac CT
 - Added the **IVUS** View for enhanced vessel visualization.
 - Added **Color Maps**, accessible through the **Show/Hide** menu.
- Vascular CT
 - Added Multi-Batch Support to improve workflow efficiency.
 - Added Spectral Data Support, for advanced spectral CT workflows.

MR applications

LoBI

- Added the functionality "Add to Running Application".

CT applications

CT Brain Perfusion (CT BP)

- Added calculation and visualization of Mismatch Ratio and Mismatch Volume.

Lung Nodule Assessment (LNA)

- Added Guideline-Based Calculations and Findings directly in Scene 1 and Scene 2.
- Relevant functionality previously available only in Scene 3 has been relocated to streamline workflow.
- Supported guidelines include:
 - Lung-RADS
 - Lung-RADS 1.1
 - Lung-RADS 2022
 - Fleischner
 - Risk Calculator
 - Additional supported guideline frameworks

Spectral applications

Spectral CT Viewer

Added a dedicated Uric Acid Workflow, optimized for Gout Disease assessment.

Multimodality applications

Multimodality Tumor Tracking (MMTT)

- Added support for Single Diameter Propagation:
 - You can propagate measurements from baseline to follow-up time points.
 - You can also propagate measurements from any follow-up time point, both forward and backward.
- Added support for propagation of Non-Measurable (Non-Target) Lesions:
 - Available from follow-up time points, with forward and backward propagation support.

Third-party applications

- Added integration with the ADAS 3D application, enabling streamlined access to external analysis workflows.
- Updated to NeuroQuant V5.3

VuePACS integration

This release introduces a tightly integrated, end-to-end workflow for Cardiovascular Suite embedded within the VuePACS Desktop Diagnostic Viewer. This enables direct access to AV Viewer, including Spectral review capabilities, Cardiac CT, Cardiac MR, and AV Vascular applications from within the VuePACS environment. Key capabilities include:

- Single Sign-On (SSO) support leveraging VuePACS authentication capabilities.
- Streamlined user experience without requiring separate application login or launch.
- Embedded launch of the Cardiovascular Suite applications from within the VuePACS viewer.
- Integrated reporting workflow between Cardiovascular Suite and VuePACS.
 - Limitation: XML reporting is currently supported only for Cardiac MR workflows.

4 Resolved issues

V16.0.3

AV Viewer

Spectral

Issue	Status
When you launched SBI data into AV Viewer, then right-clicked on Magic Glass , then selected Capture image as , the burnt label would not appear on the captured image.	Fixed. The Capture image as option is no longer available on the right-click shortcut menu.

300017891311_A/881 * 2026-08-11

5 Known issues

Note: The known issues detailed in this section that do not have workarounds (N/A) do not pose a safety risk and do not introduce any new potential security or safety hazard concerns.

Known issues from past releases that are resolved this release are moved to the "Resolved issues" section. Known issues from past releases that were already resolved before this release are deleted altogether. Due to this, the "Known issues" section changes from release to release, and is not a static record of past defects.

V16.0.3

System

Issue	Workaround
AVW does not currently support LDAPs encryption.	N/A.

Vue PACS

Known issues and Limitation in Vue PACS integration

- Selecting multiple series to launch AV Cardiovascular Suite is not supported from the VUE side. The study selector doesn't exist in this integration
 - The workaround is to select the series through the AV Cardiovascular Suite, as when the entire study is in AVW local storage, there is no fetch on-the-fly.
- XML report integration is only available for the application Cardiac MR.
- A graceful exit, where a prompt asks you to save results, only occurs when you close applications within the AV Cardiovascular Suite.
 - Otherwise, save results in the AV Cardiovascular Suite before closing the study from PACS. Any exit from Vue PACS does not trigger saving results within CV suite.
 - Auto-saving results is available only for the application Cardiac MR.
- Supported deployments:
 - HX/EX is supported.
 - Concerto and IX are out of scope.

Multimodality applications

Multimodality Tumor Tracking (MMTT)

Issue	Workaround
When you load two studies to the application, then mark a non-target lesion on a baseline study, and then on the lesion properties you mark it as Lesion disappeared , the lesion propagates to the next time point as a Lesion disappeared with a checkmark on the non-target list. When you then select the lesion properties, and deselect the checkbox Lesion disappeared , then click OK , and then either right-click the lesion or check the floating toolbar, although the " lesion disappeared " checkbox was unchecked, the option to propagate to the next time point was not available as expected.	Use the "unmatch" function for the lesions, and then propagate.

AV Cardiovascular Suite and AV Viewer

Issue	Workaround
When you zoom and pan on cMPR layouts and save these changes as a bookmark, then reload the case and activate the bookmark, the zoom and pan changes you saved do not restore.	N/A.
When you launch one study to AV Viewer, and then select Add to running to add a second study, the Label Data dialog does not appear on the second study, so there is no option to manually label it.	Exit and then re-launch the added study.
When you load any case in AV Viewer, then open AV Vascular and show both applications side-by-side by changing the display layouts to 1x2, then go to AV Viewer and change the view mode to Curve, then select an existing curve from the list, then go to User Preferences and open the AV Vascular Application, then select the vessel naming option and rename the selected vessel from the list, the AV Viewer application does not refresh the list and the old vessel name remains.	Exit curve mode and then re-enter it.

AV Vascular applications

Issue	Workaround
When you launch the application with a CT head study, then go to Volume and tissue editing, then perform free/smart sculpting in the second VR viewport (without bones VR), the vessels render in the first VR viewport (bones VR) as expected. However, when you perform free/smart sculpting in the first VR viewport (bones VR), the vessels are not rendered in the second VR viewport (without bones VR).	For Head & Neck studies, do not use free/smart sculpting when Show main vessels is enabled.
When you navigate to Volume and Tissue Editing, then select Inject Dye and proceed to inject dye into tissue, then from the upper toolbar select one of the general navigation tools, e.g., Pan, Zoom, or Windowing, and then return to the dye injection, the injection action is canceled and any unsaved injected tissue is lost. There is no option to return to the active injection.	Before using the general navigation tools, complete the Inject Dye action by selecting Confirm .
Adding private presets to a public Multi-batch group may cause a crash.	Avoid adding private presets to a public group, and vice versa.
From the Vessel list , when you select a vessel tree, then scroll on the active viewport, then select an inactive viewport, the inactive viewport becomes active as expected, but neither the image nor the contours or other annotations render.	Scrolling recovers the blank viewport, and the image renders. Maximize the viewport.
Sometimes launching a Head study fails.	N/A.
In the AVA application, when you launch a study via PACS Plugin, and then you try to perform DICOM export operations, the PACS remote node is not enabled and locked in the Save dialogs as expected.	In the Preferences, for the Save option, mark the remote node as the default save location.

AV Cardiac CT applications

Issue	Workaround
When you load a study without ECG/contrast info via PACS Plugin, then correct the labels and save your corrections, results generate and CCA launches. Upon exiting and relaunching the study, the saved labels are not retained, and the CCA launch fails.	Sometimes when launching using PACS integration, data is not identified as multi-phase cardiac, and then the Multiphase hanging does not working as expected. Open the Label dialog and update the ECG and contrast attributes for the cardiac phases.

Issue	Workaround
Cardiac CT multiphase studies with no available ECG/ Contrast DICOM info, which were originally scanned on a GE Revolution Evo, fail to trigger Cardiac CT automatic preprocessing.	Use the Preferences > Processing to Add a manual preprocessing rule. Then copy a study from PACS, and CCA launches with the manual preprocessing results.

AV Viewer

Spectral

Issue	Workaround
For a Philips Spectral Cardiac CT study containing both SBI (Spectral-Based Imaging) and non-SBI cardiac phase series, AV Viewer fails to automatically generate the Combined Multiphase Series when all phase series are selected.	To successfully generate the Combined Multiphase Series, manually select part of the series by excluding SBI-derived series, and selecting only the non-SBI phase series.
In Spectral fusion, measurements and captures do not get saved to the Findings Dashboard.	Take a screen capture for the measurement, using the Spacebar or the right-click shortcut menu.

Issue	Workaround
The Cine-activated Hanging Protocol does not load all Short Axis images, and instead shows blank viewports.	After upgrade or first use, sometimes not all the Cardiac MR orientations are automatically identified. Use the Labeling dialog as needed to manually define the orientation of these series. This is required only once.
When you load an SBI series to AV Viewer, then from QSS you select the result series Iodine with no water, then from 3D Mode you generate a batch series, then add that series to a running application, then load the batch series and press F to send the image to Film, AV Viewer crashes.	For any 3D batch created for Non-HU data types, for Cardiac multiphase data, instead of using the keyboard shortcut to Film, use the Quick Film option from the Patient Directory .
Rarely, a warning appears stating that the AV Viewer cannot be launched when more than 10 users are connected to the server, even when there are less than 10 active users.	N/A. This is a licensing issue. The Philips Service Engineer can look at all the available licenses on the Concerto CMT page. If Zero Footprint Viewer is deactivated during license generation, the AV Viewer option does not appear on the CMT license page. However, the AV Viewer and AV Applications can still launch.

Issue	Workaround
When prior studies are available on your local folder, and then you log into AVW and load a current study that has prior into AV Viewer, and the prior study is available on the study timeline, you select a 1X2 layout and then open the current and prior study side-by-side, and then right-click the image viewport, the batch controls are missing for the prior study, including Advanced batch, Quick batch, etc.	You cannot generate a batch on a study that was not originally loaded as the "Main" study. Typically this means that if you load the current study, you cannot generate a batch on its priors. Exit and select the study you need to use for batch generation, then reload. Now this study is the "Main" study and you can generate a batch on it
When you collapse the left panel of the Viewer and the Tissue list is open, right-clicking a selected tissue to open the shortcut menu hides the Tissue list.	Expand the left panel of the Viewer.
When you switch to Slab mode, create a manual curve, then switch to Curve mode, then change the inner layout to 2 CMPR, the manual curve always loads on Edit mode and causes the CMPR layout to appear blank.	If no editing is required, exit Curve mode editing by selectig the Edit curve button from the toolbar.
When you launch a study through the PACS Plugin, open the Series Selector dialog, and select a series, you cannot select Add Data to Running Applications .	N/A. When you launch through the PACS Plugin, the Add Data to Running Applications feature is not available.
The contextual toolbar is lost when you enter into an Advanced batch, close the task, and then from the leftmost panel you select the either the Cancel or Generate button.	Click again on the selected viewport, the floating toolbar appears.
When you launch a Cardiac CT Study with results, select Curve mode, then select one of the curves from the Curve list, then save a bookmark, then switch to another view mode, and then from the Findings Dashboard you activate the bookmark, a blank screen appears.	For curves that were created in the Viewer, you can activate the curve from the Findings Dashboard. For an application results curve, select the curve in the Curve mode.
When you use a PACS integration workflow, you cannot add a study to a running application if it was not pushed first to the AVW database.	Copy the study to AVW before adding it.
AV Viewer does not launch and a blank screen appears when you try to launch an NM Modality study from the Patient Directory's Images tab.	Load the study from the Series tab of the Patient Directory .


Service & IT

Issue	Workaround
When you log into the AVW server as a remote admin user: - You can access the D drive on the server (D: Portalspms/Data) - From the SQL server management studio, you can access the patient database.	N/A.
The installation process of 3mensio generates two registry entries. When the standalone installation is executed, there is only one entry visible in the Add or Remove Programs feature. However, when installation is initiated through the bootstrapper, the entry in Add or Remove Programs is enabled, which overwrites the previous setting, thereby displaying two entries on the program list.	N/A. There is no impact to the end user, and the application's launch is not affected.
Sometimes, when data exports to VuePACS fail due to connectivity issues, the job details do not appear when you select the session in Temp > VuePACS	Failed job details are stored in PiiJobServiceDB .

www.philips.com/healthcare
healthcare@philips.com



Philips Medical Systems Nederland B.V.
Veenpluis 6
5684 PC Best
The Netherlands

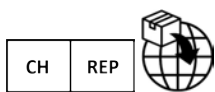

Philips AG, Seestrasse 87, CH-8810 Horgen

CE 0344



© 2026 Koninklijke Philips N.V.

All rights reserved. Reproduction or transmission in whole or in part, in any form or by any means, electronic, mechanical or otherwise, is prohibited without the prior written consent of the copyright owner.

Copyrights and all other proprietary rights in any software and related documentation ("Software") made available to you rest exclusively with Philips or its licensors. No title or ownership in the Software is conferred to you. Use of the Software is subject to the end user license conditions as are available on request.

Printed in The Netherlands

300017891311_A/881 * 2026-08-11 - en-US